
**ultraviolet shoe sanitizer  
instruction manual**

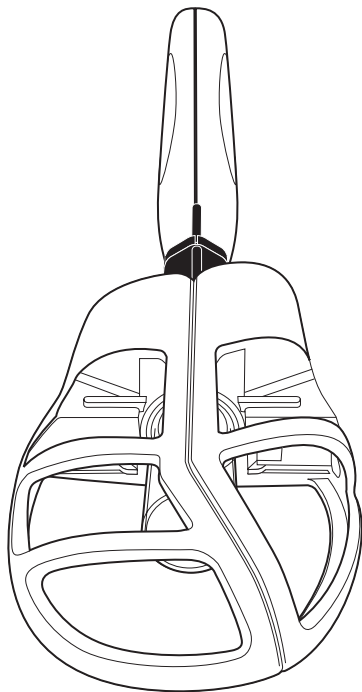**English  
Français  
Deutsch  
Español**

# IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following.

## READ ALL INSTRUCTIONS BEFORE USING

**DANGER** - To reduce the risk of electrocution:

- Always unplug the SteriShoe® shoe sanitizer immediately after using.
- Do not use while bathing.
- Do not place or store the SteriShoe® shoe sanitizer where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.
- Do not reach for a SteriShoe® shoe sanitizer that has fallen into water. Unplug immediately.

**WARNING** - To reduce the risk of burns, electrocution, fire, or injury to persons.

- The SteriShoe® shoe sanitizer should never be left unattended when plugged in.
- Close supervision is necessary when the SteriShoe® shoe sanitizer is being used by, on, or near children or invalids.
- Use the SteriShoe® shoe sanitizer only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.

- Never operate the SteriShoe® shoe sanitizer if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
- Keep the cord away from heated surfaces.
- Never drop or insert any objects into any opening.
- Do not look at the lamp when it is on.
- Do not use the SteriShoe® shoe sanitizer outside a shoe.
- Treat the SteriShoe® shoe sanitizer carefully. You are placing a fragile lamp into a shoe.
- If the SteriShoe® shoe sanitizer does not turn off after approximately 45 minutes, the timer may not be functioning properly. Turn the unit off, disconnect the power, and please call 1-866-686-Shoe(7463).
- If the on/off light is flashing, inspect the lamp to be sure that it is not broken. If the lamp is broken, **DO NOT** put your foot into the shoe until you have removed all of the pieces of the broken lamp.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.

## SAVE THESE INSTRUCTIONS

### Specifications:

SteriShoe® ultraviolet shoe sanitizer, Model: STSH-002, 20VDC, 300 mA

Power Adaptor: Input: 100 -240 V, 50/60 Hz, 0.68A Max., Class V; Output: 20VDC, 300 mA

UV lamp: 10VDC, 3W with a peak frequency of 253.7nm.

**SAA120280**

### Correct disposal of this product

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources.

To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

## Thank you for becoming a SteriShoe® customer.

You have taken an important step towards keeping your feet healthy.

### The SteriShoe® sanitizer is:

- Scientifically proven to destroy microorganisms.\*
- Effective on the microorganisms that cause toenail fungus, athlete's foot, and offensive foot odor.
- Non-toxic. Uses no harsh chemicals, sprays or drugs to kill bacteria, fungi, and germs.
- Accepted by the American Podiatric Medical Association.
- Very easy, safe and convenient to use. Just plug it in and push a button!

Please read through this manual to learn how to use your SteriShoe® shoe sanitizer safely and effectively. Let us know how you are enjoying your SteriShoe® sanitizer by e-mailing us at [feedback@SteriShoe.com](mailto:feedback@SteriShoe.com). We love to hear from our customers.

Sincerely,  
Adam Ullman  
President & CEO

\*Independent testing results available at [www.SteriShoe.com](http://www.SteriShoe.com)

## Use your SteriShoe® sanitizer daily

The combination of sweat and heat inside your shoes makes a perfect breeding ground for the bacteria and fungi that cause toenail fungus, athlete's foot, and foul odor. And when foot ailments exist, shoes can be a common source of re-infection.

Therefore, to fully minimize the likelihood of a new infection or re-infection, germs must be eradicated at the risk source daily.

We recommend that you change the UV lamps in your SteriShoe® every two years to guarantee maximum effectiveness killing harmful microbes in your shoes. We will make it easy for you to remember.

**Register your SteriShoe® shoe sanitizer at [www.sterishoe.com/regist.html](http://www.sterishoe.com/regist.html). We will send you an email when it is time to replace the UV lamps.**


## How the SteriShoe® sanitizer works:

The innovative SteriShoe® shoe sanitizer kills the germs that reside inside shoes with a special germicidal ultraviolet light (UVC).

With the first treatment, a SteriShoe® sanitizer destroys up to **99.9%** of the microorganisms in shoes, including the bacteria and fungi that cause toenail fungus (onychomycosis), athlete's foot (tinea pedis) and offensive foot odor.

The SteriShoe® sanitizer is also an excellent choice for those with diabetes or generally concerned about their feet's hygiene.

## Contents
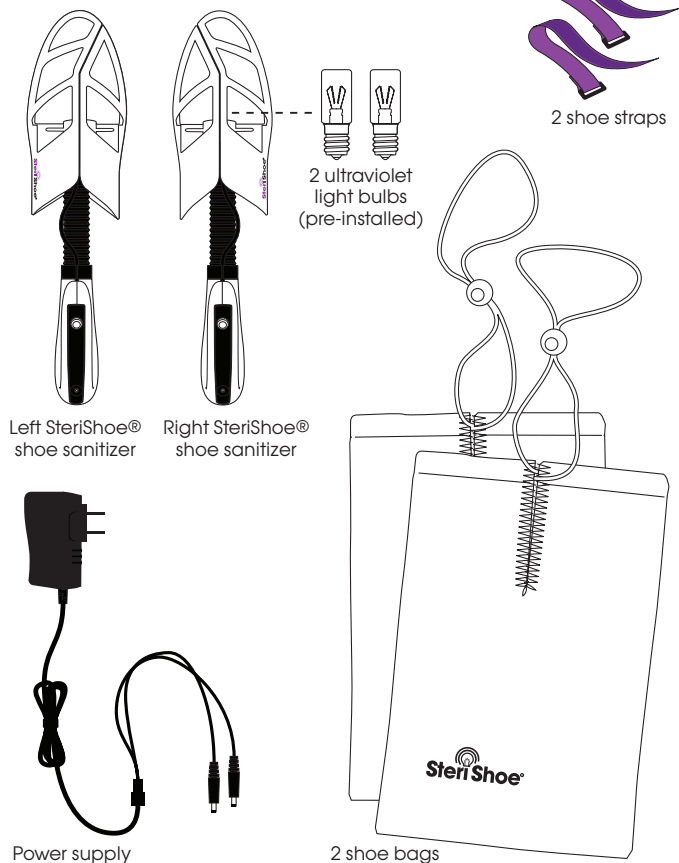
## Using your SteriShoe® device

Your SteriShoe® sanitizer comes with important safety features which prevent exposure to the UVC light. The following is a quick reference on how to safely operate your SteriShoe® ultraviolet shoe sanitizer.

The on/off button functions as a status indicator.

Status states:


no light =  
unit is off

solid light =  
unit is treating a shoe

blinking light =  
an error has occurred

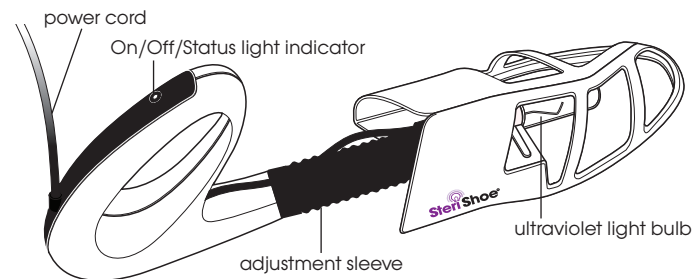
## Use the SteriShoe® shoe sanitizer only as directed.

To prevent errors when using the SteriShoe® shoe sanitizer, please follow these directions:


### Installing or replacing a lamp

Using a soft dry cloth, hold the front of the UV lamp and screw it into the socket in the SteriShoe® ultraviolet shoe sanitizer. Repeat this step for the second SteriShoe® shoe sanitizer.


### General Use Instructions

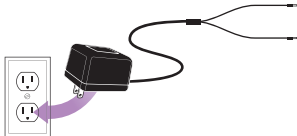
1. Plug the power supply into a standard wall outlet (AC 100-240V).

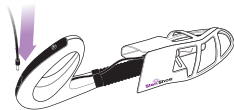
2. Plug one power cord into each SteriShoe® shoe sanitizer.

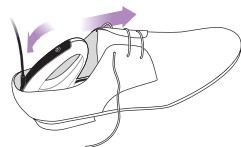3. Insert the right SteriShoe® shoe sanitizer into a right shoe and the left SteriShoe® shoe sanitizer into a left shoe.

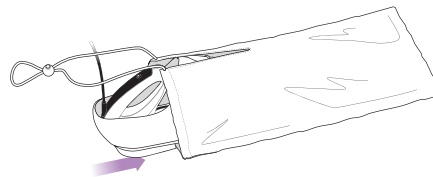4. Place each SteriShoe® shoe sanitizer and shoe into the provided SteriShoe® branded shoe bag.

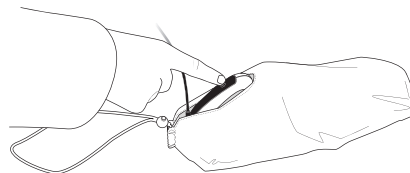5. Close each shoe bag drawstring.

6. Press each on/off button to begin treating your shoes.

7. If the on/off button is glowing solid purple, the SteriShoe® shoe sanitizer is treating the shoe and will turn off after approximately 45 minutes. **DO NOT** remove the shoe bag while the SteriShoe® shoe sanitizer is treating a shoe.

8. Once the treatment process is complete, the on/off light will no longer be illuminated and it will be safe to remove the SteriShoe® shoe sanitizer from a shoe.

If you need to remove a SteriShoe® shoe sanitizer from a shoe when it is treating a shoe, press the on/off button to turn it off. If the on/off button is not solidly illuminated, then it is safe to remove the SteriShoe® shoe sanitizer as it will not interrupt a treatment.

**Note:** If the on/off button is flashing, one of the safety features was probably activated. See the troubleshooting section (on page 13) for more information.


## Using your SteriShoe® device

The SteriShoe® ultraviolet shoe sanitizer will not work if:

- You use it in a shoe without heel support. The SteriShoe® shoe sanitizer must be compressed to operate. (See instructions for using the SteriShoe® shoe strap on the opposite page.)
- It detects ambient light. The provided shoe bag must be used to darken the shoe environment.

Just like you give your feet a break, give your SteriShoe® shoe sanitizer a break too. Allow your SteriShoe® shoe sanitizer to rest for at least fifteen minutes between uses.

## Using your SteriShoe® shoe strap

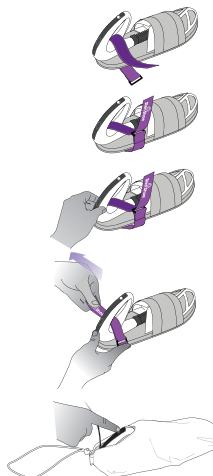1. Insert opened strap through the SteriShoe® shoe sanitizer handle (with the Velcro material facing outward)
2. Insert strap through the d-ring on the end of the SteriShoe® strap.
3. Push the heel of the SteriShoe® shoe sanitizer forward.
4. While pushing the heel forward, tighten the strap and secure it with the Velcro.
5. Plug in your shoe sanitizers, insert them into the SteriShoe® shoe bags and turn them on.

**DO NOT** allow the shoe strap to go between the front halves of your SteriShoe® sanitizer. For use with shoes where this will occur, reverse the shoe sanitizers (i.e. use the right in the left shoe and the left in the right shoe) to prevent this issue.

## Troubleshooting

There are safety features in the SteriShoe® shoe sanitizer to prevent a user from being exposed to the UV light. If the on/off button is blinking, one of the safety features has likely been activated. The error chart is listed below. If you have any questions, please call 1-866-686-7463, ext 2 for technical assistance.

### Error and Remedy

The ambient light sensor has detected light. Be sure that the SteriShoe® shoe sanitizer is inside a shoe. If the shoe is a sandal or open shoe, be sure to place the SteriShoe® shoe sanitizer and shoe inside the provided shoe bags. If the error is still present, try pressing down on the tongue of the shoe when turning the SteriShoe® shoe sanitizer on or move to a darker room.

The SteriShoe® shoe sanitizer is either not in a shoe or the shoe is too large to compress the safety switch. Try using the SteriShoe® shoe sanitizer in a smaller shoe or

place something, such as a pair of socks, between the heel of the shoe and the SteriShoe® shoe sanitizer. If you are using the SteriShoe® shoe sanitizer with a shoe that does not have heel support, use the SteriShoe® shoe strap to keep the SteriShoe® sanitizer compressed.

The lamp is either damaged or not screwed in completely. If you have this error, unplug the power and switch the lamps between the left and right shoe sanitizers.

- If both units now work (the on/off light remains illuminated), the lamp may have not been screwed in completely.
- If the problem remains with the lamp, e.g. the right shoe sanitizer was not working and after switching the lamps the left shoe sanitizer is not working, then the lamp needs to be replaced.

**Note:** After turning the shoe sanitizer on, it will take approximately 10 seconds for this error indicator to occur. This is the time needed to ignite the UVC lamp.

## Important Information

DO NOT look at the light from your SteriShoe® ultraviolet shoe sanitizer when it is turned on; it emits UV light which can be harmful to eyes and skin. You

can tell when the SteriShoe® shoe sanitizer light is on because the on/off button will be glowing solid purple.

If the on/off button is flashing, an error has occurred. The SteriShoe® shoe sanitizer has built-in safety features to reduce the likelihood of accidental exposure to the UV lamp. An error is most likely due to the activation of a safety feature. See the troubleshooting section for more information.

**LAMP CONTAINS MERCURY (Hg).** The lamp contains a tiny metal ball which you may hear inside the lamp. This is normal; the lamp is not broken. Since mercury is toxic please dispose of the lamp in accordance with local disposal laws. For more information, see: [www.lamprecycle.org/res-statelaws.shtml](http://www.lamprecycle.org/res-statelaws.shtml).

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

**Warning:** Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

### Disclaimers

- The SteriShoe® shoe sanitizer destroys up to 99.9% of the targeted microorganisms. Targeted microorganisms are those that come in direct unobstructed contact with the UV light.
- The SteriShoe® shoe sanitizer does not treat medical conditions; it is intended to provide a clean and healthy environment for feet inside shoes.
- The UV light in the SteriShoe® may discolor the interior of the shoe.
- UV light may accelerate the degradation of certain materials.

## Warranty Information

### Limited One Year Warranty

**Warranty coverage:** Shoe Care Innovations, Inc. ("SCI") warrants, to the original purchaser, that this product will be free of defects in materials or workmanship for a period of one (1) year\* from date of purchase, with the exceptions stated below.

\*The warranty term for lamps is 60 days from the date of purchase.

**What is not covered:** This warranty does not cover damage caused by accident, tampering, improper care, misuse, uses for which the product was not designed, or other causes not arising from defects in material or workmanship.

**Warranty performance:** During the above warranty period, SCI will repair or replace any part of the product deemed defective.

Contact SCI at 1-866-686-Shoe(7463) to obtain warranty service or go to [www.SteriShoe.com](http://www.SteriShoe.com) for more information.

**Warranty disclaimers:** Any implied warranties arising out of this sale, including but not limited to the implied warranties of description, merchantability and fitness for a particular purpose, are limited in duration to the stated warranty coverage. In no event shall SCI be liable for loss of use of this product or for any indirect, special, incidental or consequential damages, or costs, or expenses incurred by the consumer or any other user of this product, whether due to a breach of contract, negligence, strict liability in tort or otherwise. SCI shall have no liability for any personal injury, property damage or any special, incidental, contingent or consequential damage of any kind resulting from the use of this product. Some states do not allow limitations on how long an implied warranty lasts or the exclusion or limitation of consequential or incidental damages, so the above limitation may not apply to you.

**Legal remedies:** This warranty gives you specific legal rights and you may also have other rights that vary from state to state.

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Copyright © 2012 Shoe Care Innovations, Inc.
SteriShoe, the SteriShoe logo, and the Shoe Care Innovations, Inc. logo are registered trademarks of Shoe Care Innovations, Inc. in the U.S. and other jurisdictions.